
# ONBOARDING GUIDE

---

## Withings Body Pro

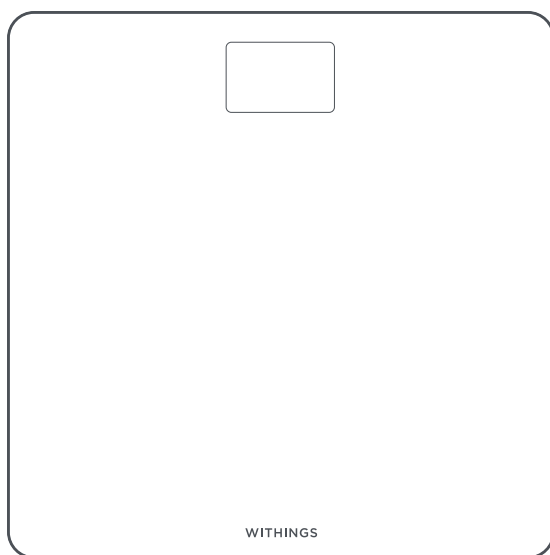Guide d'installation | Guía de instalación |  
Installationsanleitung

EN | FR | ES | DE

**WITHINGS**

---

THIS DEVICE HAS BEEN  
PROVIDED TO YOU BY  
YOUR HEALTH MONITORING  
PROGRAM.

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
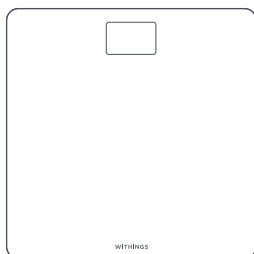
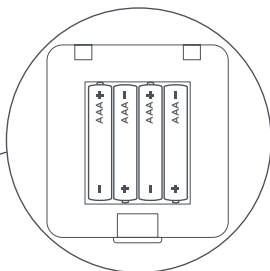
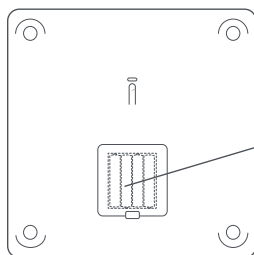
Alkaline cells  
already installed



Regulatory information leaflet

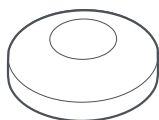Carpet feet (optional installation)

## OPTIONAL : Installation of carpet feet

Use only on soft surfaces as carpet feet can impact accuracy if used on hard floors.
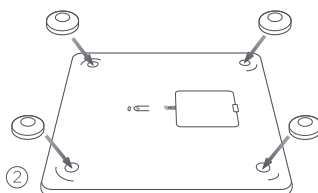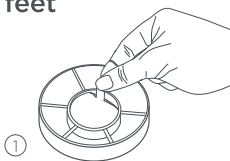
# Installation



Body Pro scale is ready to install right out of the box.

Setup can take a few minutes and only happens the first time you use your scale.

Place the scale on a hard, flat surface. If you have carpet, please proceed to carpet feet installation.

## ① Wake up the scale

Pull the tab underneath the scale that is found inserted in the battery cover.

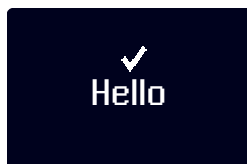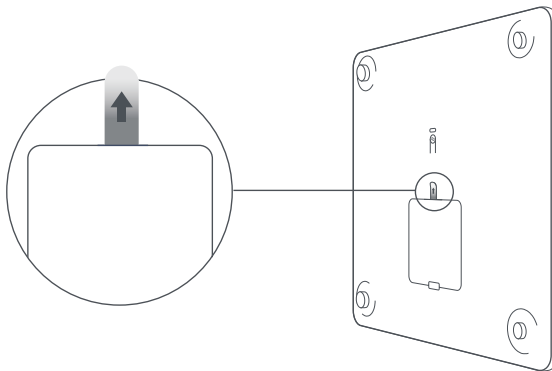
The screen will display **“Hello”** when it wakes up.

## ② Wait for the scale to set up

Body Pro will set up automatically in three steps. You can follow the progress on the scale's screen.

The installation process can take up to 15 minutes. Please do not move the scale during setup.


The scale will display **“Ready”** on its screen once setup is complete.
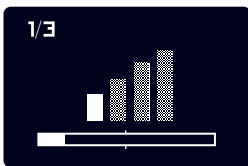
## ③ Take your first measurement


Step on the scale to weigh yourself.

At your first weigh-in, the scale will display **“Setup successful.”**
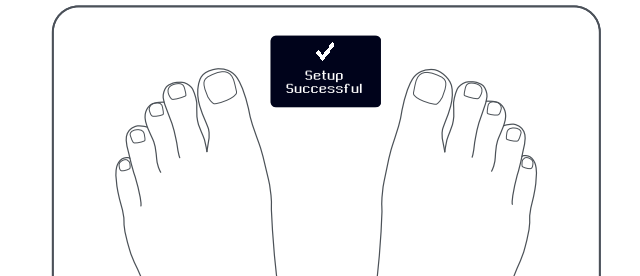
## Best practices

- If you move your scale, wait for 2 seconds before stepping on it.
- Step on the scale with clean, bare feet.  
Wear as little clothing as possible (or no clothes at all).
- Do not bend your knees when weighing yourself.
- To ensure a highly accurate measurement, arrows may appear on your screen to help you distribute your weight evenly on the scale's surface. If one or more arrows appear on any corners of the scale's display, simply shift your weight in the opposite direction of the arrows until they disappear.

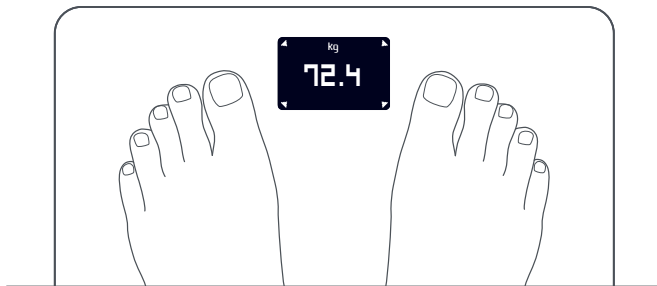
**Position Control**

- Try to weigh yourself at the same time each day. For example, make it a habit to step on in the morning before you eat or drink anything.
- When you step off the scale, it will power off automatically.

You can change the units of measure on the Body scale (from kg to lb or st lb). Follow this procedure to do so:

- Give a quick press on the button at the back of the scale. This will show the three options for the units of measure (kg, lb, and st lb) on the screen.

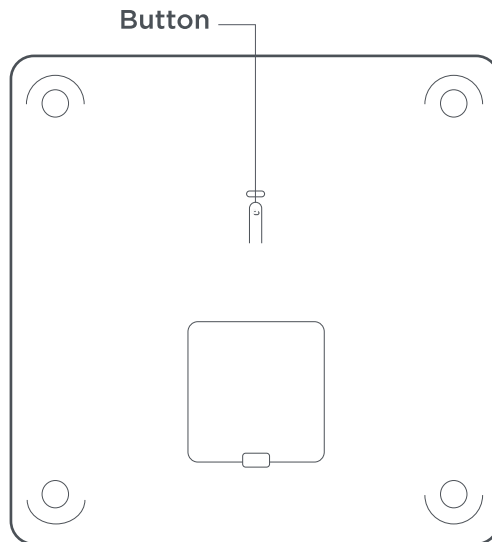- Press the button again to move the arrow on the screen to the desired unit. When you are on the desired unit, confirm with a long press. The units of measure have now been changed.

# Get to know Body Pro


Along with weight, see your Body Mass Index and the local weather forecast to help plan your day while weighing in.

Weight (kg, lb, st lb)

BMI

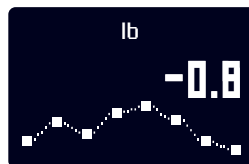Weight trend (kg, lb, st lb)

Weather (morning / afternoon,  
°C or °F depending on location)

## During the installation

| PROBLEM   | REMEDY  |
|---|---|
|  | <p>Please try installing Body Pro again by stepping on the scale or giving a quick press on the button at the back of the scale.</p>  |
|  | <p>Try changing the location of your scale by moving it close to a window where the network coverage is better, and then step on the scale or give a quick press on the button at the back of the scale to start the installation process again. If the problem persists, the issue might be due to a cellular dead zone.</p> |

## During the installation or usage

| PROBLEM                                  | REMEDY   |
|--|--|
| <p>ERROR 401</p>                         | <p>This means Body Pro was not able to upload your previous measurements and is having issues connecting to a network. Please try to change the location of your scale close to a window where the network coverage is better, and take another measurement. If the issue persists, please contact the customer service number of your health program in which you are enrolled.</p> |
| <p>ERRORS<br/>402, 403,<br/>404, 407</p> | <p>If you receive any of the error messages listed below, please contact the customer service number of your health program in which you are enrolled.</p>   |

### Need help?

Please contact our support program-support@withings.com

## Caring for your device

Withings devices are designed for long-term use. Withings will continually improve your scale even after installation. We are committed to providing you with the most accurate measurements possible now and over the long term.

To help ensure the longevity of your device, please note the following care and maintenance tips:

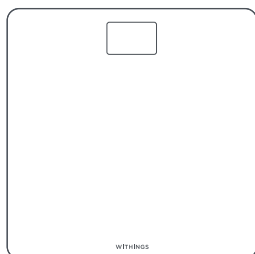
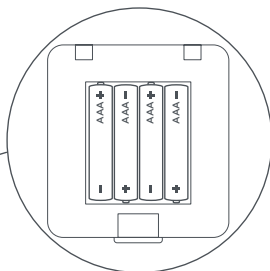
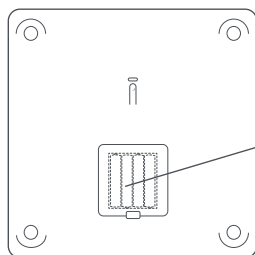
- To clean your scale, only use a cloth dampened with water. Avoid chemical cleaners, abrasive cleaners, solvents, and insect repellents that can damage plastic components and finishes.
- To prevent mishaps, do not place the scale on slippery surfaces, and do not use it when your feet or the scale is wet.
- Place your scale on a flat surface, and if necessary, make sure you are using the carpet feet to ensure accurate readings.

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© 2021 Withings. All rights reserved V2.1 - May, 2021